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Dräger Medical AG & Co. KG

For Your Safety and that of your Patients

Strictly follow the Instructions for Use

Any use of the apparatus requires full understanding and strict observation of these instructions. The apparatus may only be used for the purposes specified here.

Maintenance

The apparatus must be inspected and serviced (and a record kept) by trained service personnel after every 6000 hours of operation or at least once per year.

Repair and general overhaul of the apparatus may only be carried out by trained service personnel. We recommend that a service contract be obtained with DrägerService and that all repairs also be carried out by them.

Only authentic Dräger spare parts may be used for maintenance. Observe the chapter "Maintenance intervals".

Accessories

Only the accessories specified in the Order List may be used with the apparatus.

Liability for proper function or damage

Liability for the proper function of the apparatus is irrevocably transferred to the owner or operator to the extent that the apparatus has been improperly serviced or repaired by personnel not employed or authorized by DrägerService or if it has been used in a manner not conforming to its intended use.

Dräger cannot be held responsible for damage caused by non-compliance with the recommendations given above. The warranty and liability provisions of the terms of sale and delivery of Dräger are likewise not modified by the recommendations given above.

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Safety instructions

Not for use in areas of explosion hazard

This apparatus is neither approved nor certified for use in areas where combustible or explosive gas mixtures are likely to occur.

Safe connection with other electrical equipment

Electrical connections to equipment which are not listed in these Instructions for Use should only be made following consultations with the respective manufacturers or an expert.

Supplying life-supporting ventilators

If the compressor is used to supply air to life-support ventilators, an adequate alternative supply of medical air must be ensured in case a fault develops in the compressor!

Where the Compressor is used to supply air to life support ventilators, the ventilator should monitor the supply pressure and appropriately alarm if necessary.

Must not be used in the presence of flammable gases or anaesthetics. Fire hazard!

Inflammable liquids must not be sprayed near the compressor. Fire hazard!

Prevent pollution of room air!

The compressor draws in air from the room. Any pollutants in the air would reach the patient.

General information on electromagnetic compatibility (EMC) according to the international EMC standard IEC 60601-1-2: 2001

Medical electrical equipment needs special precautions regarding electromagnetic compatibility (EMC) and needs to be installed and put into service according to the EMC information provided in the technical documentation available from DrägerService upon request.

Portable and mobile RF communications equipment can affect medical electrical equipment.

Pins of connectors identified with the ESD warning symbol shall not be touched and not be connected unless ESD precautionary procedures are used. Such precautionary procedures may include antistatic clothing and shoes, the touch of a ground stud before and during connecting the pins or the use of electrically isolating and antistatic gloves. All staff involved in the above shall receive instruction in these procedures.

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Intended Use

Air compressor supplying compressed air for medical ventilators.

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Preparing for use

Before first use

The compressor is supplied without trolley. It can be installed on the Evita trolley or on a separate trolley.

The medical air compressor must not be started up until it has been installed on the trolley.

The compressor unit may be damaged as a result of poor air circulation if this is not done!

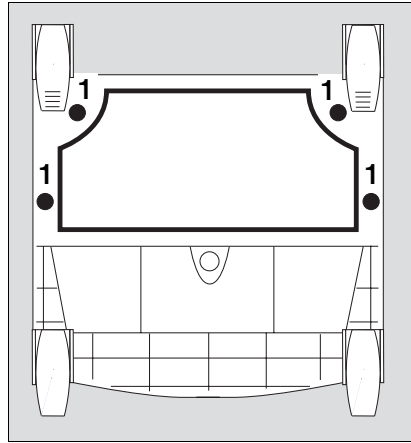
The compressor may only be installed by trained service personnel.

Tools required:


Spanner for size 5 Allen screw,
Phillips screwdriver, size 2,
Size 19 fork wrench.

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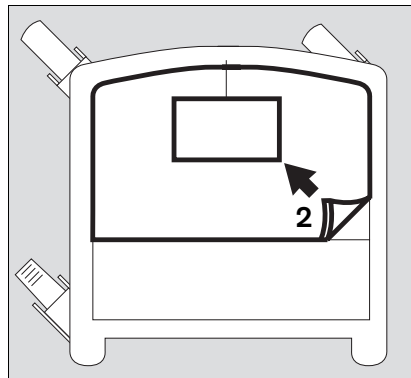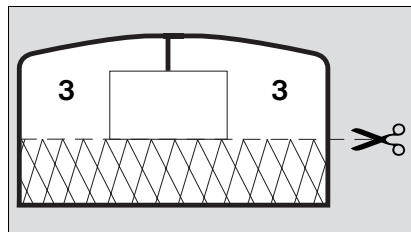
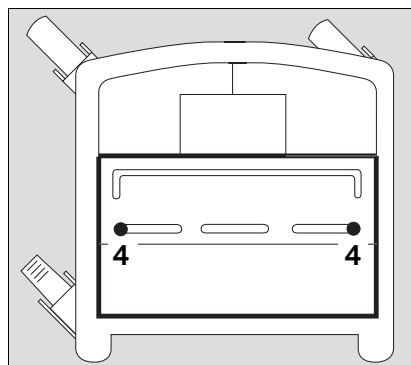
Installing the compressor on the Evita trolley

- Remove the ventilator from the Evita trolley.
- Lay the Evita trolley on the ground so that the underside is accessible.
- 1 Secure the exhaust air cover to the trolley frame with the enclosed Philips screws.
- Turn the trolley upright again.

- 2 Remove the foam mat from the trolley.

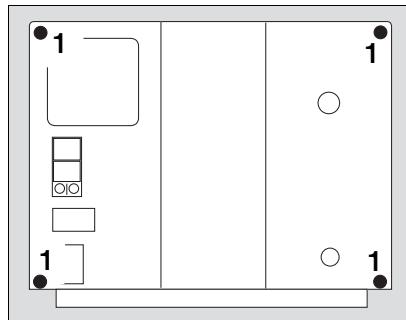
- Cut the foam mat off along the prepared line.
- 3 Place both parts in the trolley.

- 4 Undo the two screws on the battery compartment cover and remove the cover.
- Undo the two screws below the battery compartment cover.

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
On the compressor:

- 1 Undo the four Allen screws on the rear of the compressor.


Ensure that the mains power cable has been unplugged in order to avoid electric shocks!

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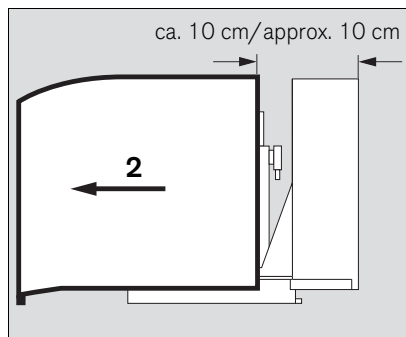
- 2 Carefully draw the shroud back approx. 10 cm so that the indicator lead is not pulled off.

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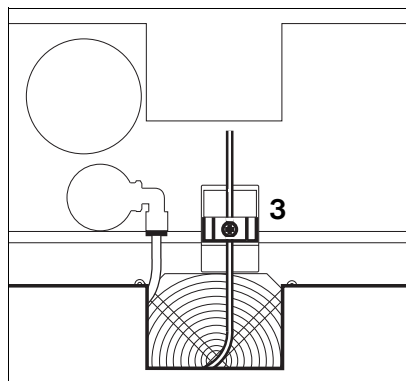
- 3 Undo the Phillips screw on the plug connector, unplug the connector and remove the shroud completely.

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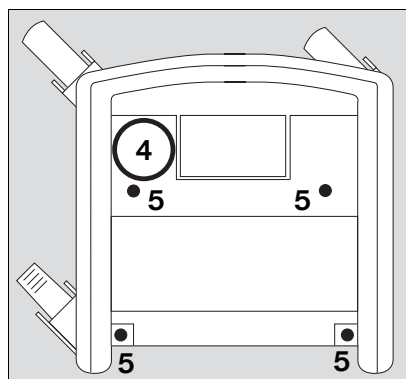
- 4 Place the compressor on the trolley with the pressure vessel alongside the column.
- 5 Tightly secure the compressor to the trolley with the four Allen screws and washers.

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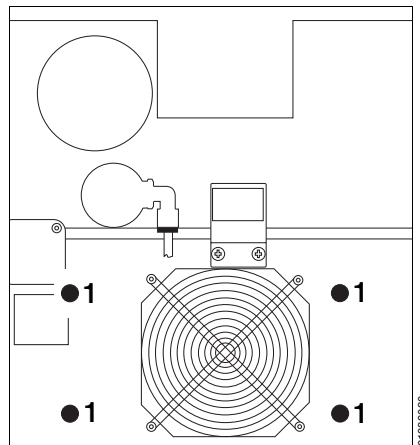- 1 Undo the four Allen screws securing the internal packing.

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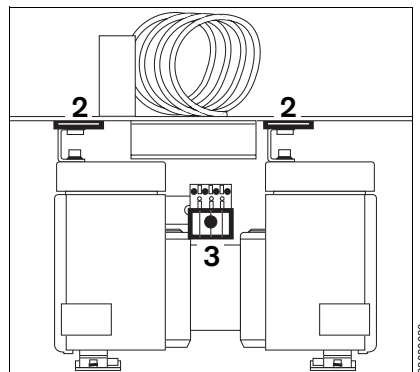- 2 Pull out the two tabs on the internal packing.
- 3 Plug in the motor connector and secure it with the Phillips screw. Keep the tabs and screws from the internal packing for future use.

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This safety mechanism must always be removed before using the compressor for the first time, otherwise major vibrations may develop in the compressor. These may damage the compressor unit!

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The internal packing must be fitted whenever the compressor is transported without its trolley, otherwise the rubber feet may be damaged.

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- Fit the shroud, plug in its connector and secure it with the Phillips screw.
- Close the shroud and screw down the four Allen screws on the rear of the unit.

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Installing the compressor on a separate trolley

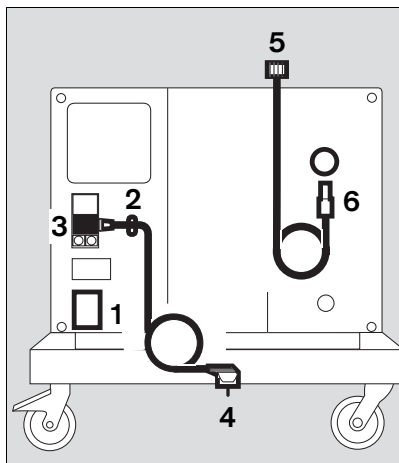
The compressor is installed in the same way as on the Evita trolley.

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Connecting the compressor

- 1 Check that the operating voltage specified on the rating plate matches the mains voltage.
- 2 Unscrew clamp.
- 3 Plug in apparatus connector, pull cable until taut, secure under clamp and
- 2 screw clamp back into place.
- 4 Plug the connector into an earthed socket; do not use an adapter.
- 5 Screw the pressure hose onto the ventilator.

**Do not connect more than one ventilator to the compressor!
It could overload the compressor.**

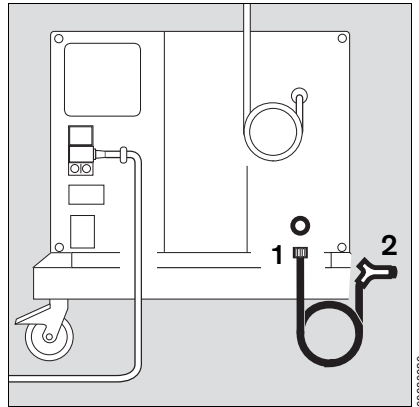
- 6 Firmly insert the connector on the pressure hose into the coupling on the compressor – until it engages.
- Place the compressor in a cool place, away from radiators and other sources of heat.
Ensure that air can circulate freely around the compressor.
- The compressor must not be operated in salty air, as this may cause the non-return valves in the compressor to corrode.

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
- 1 Screw the pressure hose to the compressor.
- 2 Plug the connector into the wall socket of the central medical air system.

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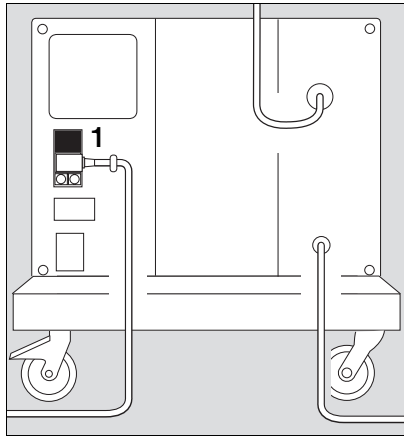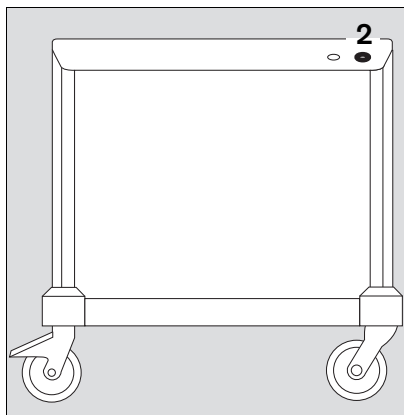
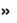
Operation

If the compressor has been stored at temperatures below 3 °C:


- Wait approx. 2 hours for the compressor to warm up.
- 1 Press main switch on rear,
- 2 the green LED »  « in the front of the shroud lights up = ON.

- Now switch on the ventilator.

Refer to "Fault - Cause - Remedy" on page 19 if any faults develop.

Standby mode

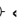For emergency supply if the central medical air system fails.

- Always leave the main switch ON = green LED »  « on front of shroud lights up.

The pressure in the central medical air system is monitored by the compressor while switched on:

If the pressure drops below 2.7–0.3 bar¹⁾, the compressor automatically supplies medical air to the ventilator to which it is connected. The compressor switches off the supply and remains on standby when the pressure in the central medical air system rises above 3.4 bar²⁾.

Switching off

After switching off the ventilator:

- Press the main switch on the rear,
- the green LED »  « goes out = OFF.

The compressor is vented automatically.

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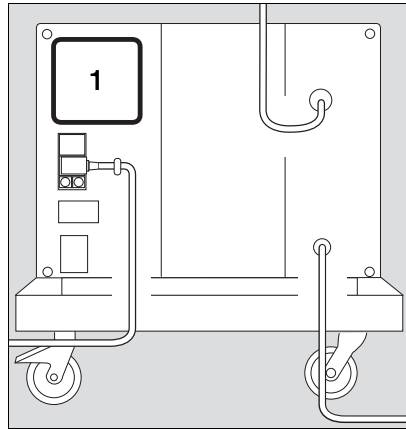
Routine maintenance

If the filter mat is soiled:

- 1 Remove the filter mat, knock out any dirt and insert the mat again.
- Any dirt on the compressor housing should be wiped off with a disposable cloth.

Wiping with disinfectant

A surface disinfectant can be used to wipe the housing. Preparations based on:

- aldehydes or
- quaternary ammonium compounds.

Preparations based on the following compounds should not be used as they may damage the housing:

- Compounds containing phenol,
- Compounds giving off halogens,
- Strong organic acids,
- Compounds giving off oxygen.
- Wipe with disinfectant e.g. Incidur® (Ecolab Deutschland GmbH)
- Note the manufacturer's instructions.

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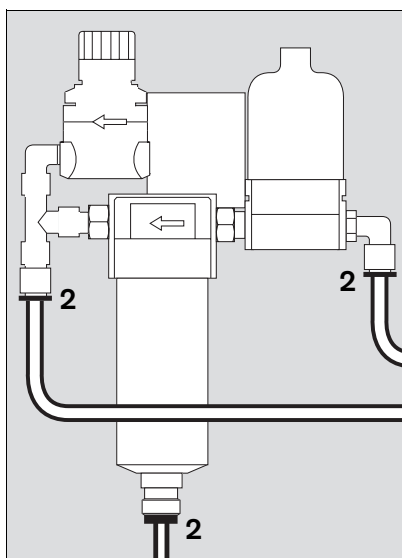
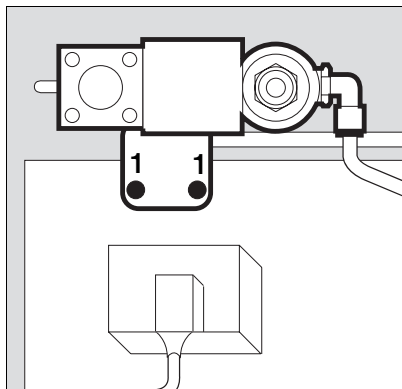
Maintenance intervals

Maintenance work may only be carried out specially trained service personnel.

at least every

6000 hours:

(see operating hours counter)

Change filters
in prefilter,
main filter and
intake filter

Carry out maintenance more frequently if necessary.

If the apparatus has not been in operation for 6000 hours:

- Carry out inspection and servicing after not more than a year.

Dispose of used filters with domestic waste.

The apparatus must always be unplugged from the mains in order to avoid electric shocks!

Removing the filter group

Remove the housing, see page 7.

- 1 Undo both Phillips screws.

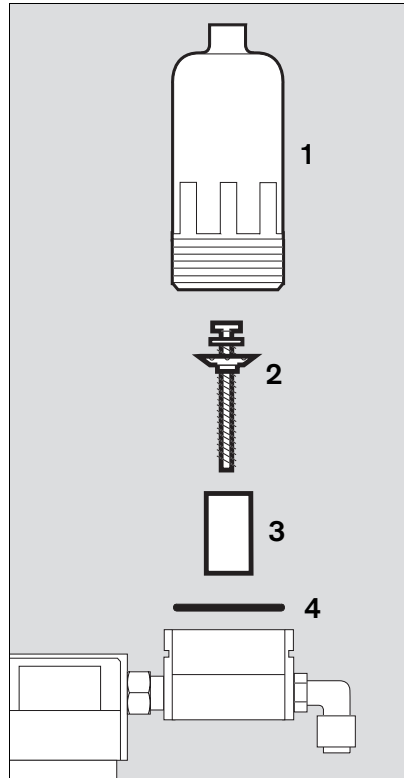
- 2 Disconnect the hoses:
Pull back and hold the ring while pulling the hose from the socket at the same time.
- Remove the filter group.

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Changing the prefilter

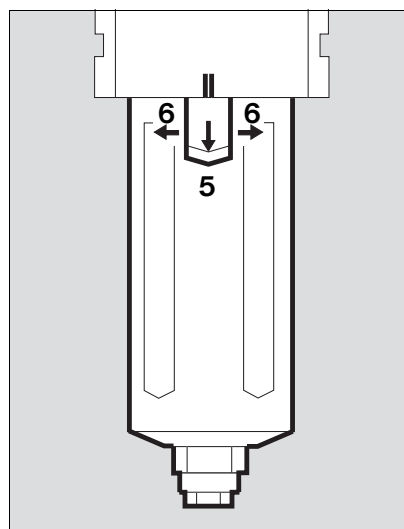
- 1 Unscrew the housing by hand.
- 2 Unscrew recessed head screw and remove screw and plastic washer.
- 3 Remove the filter sleeve and replace.
- 3 Insert the new filter sleeve.
- 2 Fasten in place with recessed head screw and plastic washer.

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
Changing the main filter


- 5 Pull and hold the catch while simultaneously
- 6 turning the housing until the marks (II) line up.
- Remove the housing.

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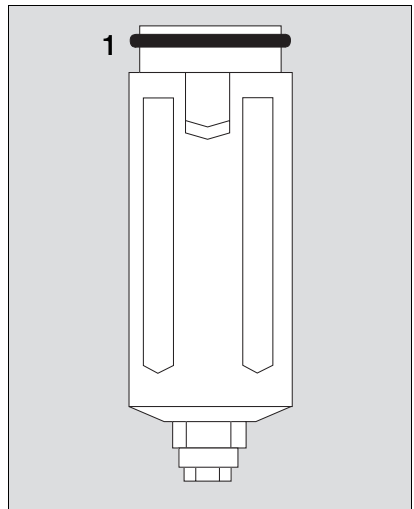- 1 Remove old O-ring from housing and fit new O-ring.

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
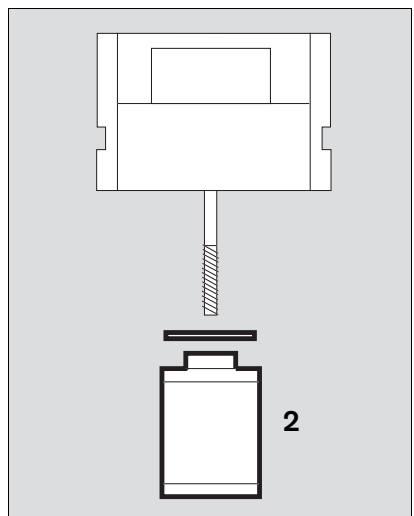
- 2 Unscrew old filter sleeve by hand and screw in new filter sleeve.
- Insert the housing so that the marks line up, then turn until the catch is heard to engage.
 - Lightly pull the housing to ensure that it has engaged.

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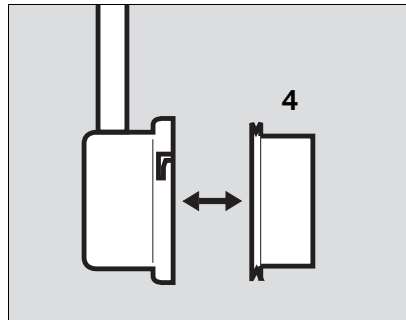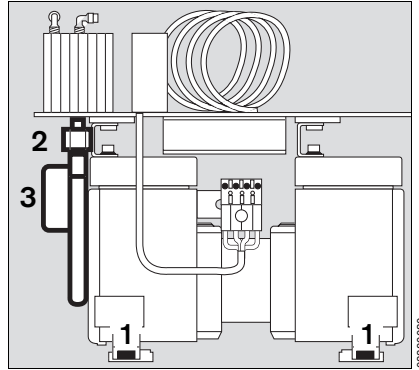Changing the intake filter

- 1 Remove the two Allen screws on the front of the assembly plates of the compressor unit and undo the two rear Allen screws two turns.
- 2 Unscrew the pressure hose with a size 19 fork wrench and hold the connector to the heat exchanger with a second spanner to prevent it working loose.
- Pull the compressor unit out until the outlet pipe at the rear of the filter housing no longer protrudes from the hole.
- 3 Turn the filter housing clockwise and remove it.
- 4 Remove the old intake filter and insert a new filter in the filter housing.
- Position the filter housing with the pipe facing to the rear in a position corresponding roughly to "11 o'clock".
- 3 Turn the filter housing anticlockwise until locked. The pipe is now horizontal.
- Pull gently to check that the filter housing is secure.
- Push the compressor unit back, guiding the pipe into the hole at the same time.
- 2 Secure the pressure hose with a size 19 fork wrench and hold the connector to the heat exchanger with a second spanner to prevent it working loose.
- 1 Secure the compressor unit with the front and rear Allen screws.

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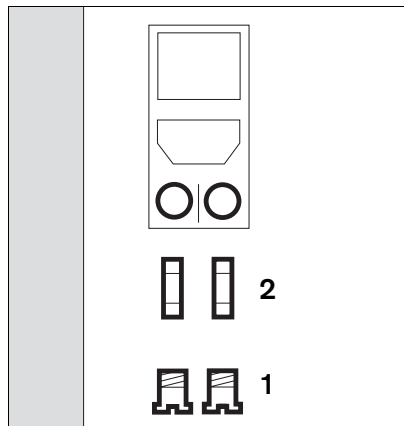
Installing the filter group

- Connect the hoses to the filter group, see page 13.
- Screw the filter group into place, see page 13.
- Fit the shroud, plug in its connector and secure it with the Phillips screw.
- Close the shroud and screw down the four Allen screws on the rear of the unit.

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Replacing the fuses

- Unplug the mains connector.
- 1 Turn the caps anticlockwise with a screwdriver and remove them.
 - 2 Remove blown fuses from their caps and insert new fuses. Only the specified fuses may be used, see "Technical data" on page 21.
 - 1 Screw the caps with the new fuses back into place.

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Disposal of the device

- at the end of its useful life.
- After contacting the competent specialised waste disposal companies, send the device to the appropriate expert disposal facility. Observe all applicable legal requirements.

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This device is subject to EU Directive 2002/96/EC (WEEE). It is not registered for use in private households, and may not be disposed of at municipal collection points for waste electrical and electronic equipment.

Dräger Medical has authorized a firm to dispose of this device in the proper manner: for more detailed information, please contact your local Dräger Medical organization.

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
¹⁾ US-Kompressor 115 V/60 Hz/Standby: 2,7±0,2 bar

Fault - Cause - Remedy

| Fault | Cause | Remedy |
|---|--|--|
| Compressor vibrates strongly | Internal packing has not been removed | Remove internal packing, page 8 |
| | Rubber feet on compressor unit defective | Call DrägerService |
| Compressor runs, but builds up too little or no pressure | Intake filter soiled | Replace intake filter, page 15 |
| | Leaks in the compressor components | Call DrägerService |
| Red LED »Temp  « lights up and an alarm sounds | High ambient temperature | Place the compressor in a cool place and ensure free air circulation |
| | Filter mat soiled | Knock dirt out of filter mat or fit a new filter mat |
| | Fan defective | Call DrägerService |
| Compressor refuses to start | Mains voltage too low or not available | Mains voltage too low or not available |
| | Fuses defective | Replace fuses, page 17 |
| | Apparatus defective | Call DrägerService |
| Standby apparatus: | | |
| Compressor starts up although pressure in central medical air supply is over 2.7-0.3 bar ¹⁾ | Pressure switch defective | Pressure switch defective |

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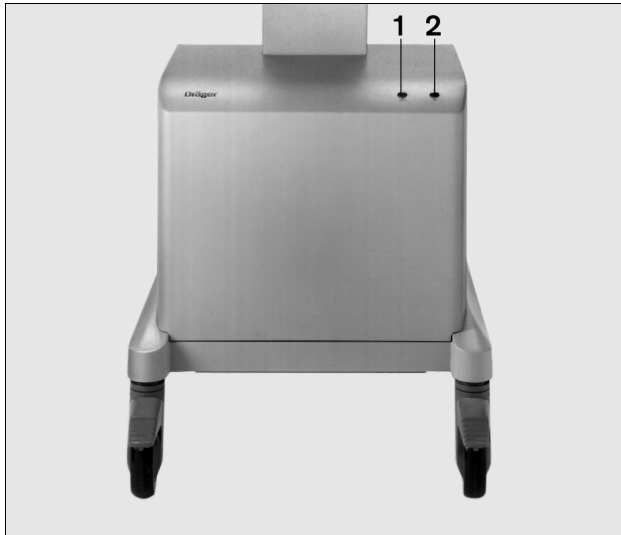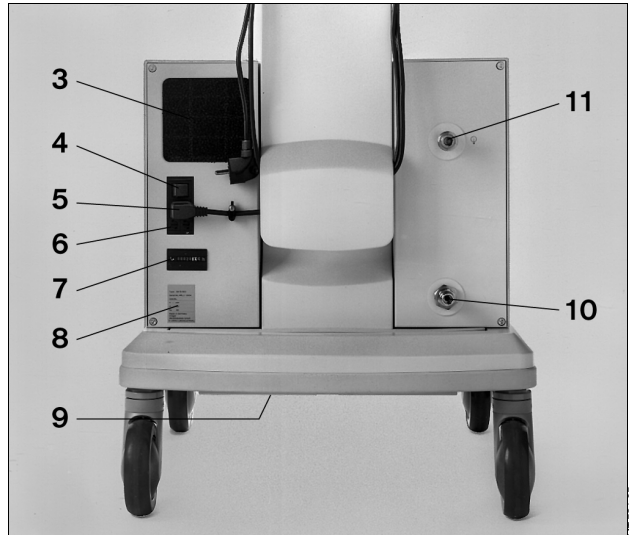
What's what

- 1 Red LED »Temp Δ« when temperature is too high
- 2 Green LED »  « indicates power on / standby
- 3 Filter mat
- 4 Main switch
- 5 Power cable
- 6 Fuses
- 7 Operating hours counter
- 8 Rating plate
- 9 Exhaust air duct
- 10 Connection for pressure hose from central medical air supply (standby apparatus only)
- 11 Coupling for pressure hose from ventilator

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Ambient conditions

Operation:

| | |
|------------------------|-----------------------------|
| Temperature | 10 to 40 °C |
| Relative humidity | 30 to 95 %, no condensation |
| Height above sea level | 0 to 4000 m |

Storage:

| | |
|------------------------|----------------------------|
| Temperature | -20 to 70 °C |
| Relative humidity | 0 to 99 %, no condensation |
| Height above sea level | 0 to 16000 m |

Performance characteristics

| | |
|---|---|
| Operating pressure | 4+0.5 bar |
| Continuous flow | Minimum 30 L/min at 3.0 bar |
| Peak flow | 180 L/min for max. 0.8 seconds at a total of max. 30 L/min |
| Dew point depression under operating pressure | 5 °C below ambient temperature at ≥ 30 L/min and max. ambient temperature 40 °C. 15 °C below ambient temperature at 15 L/min and max. ambient temperature 40 °C |
| Air quality | Dust-free and oil-free compressed air |
| Filter | $\leq 1 \text{ micron}$ |
| Compressed air outlet | Coupling with check valve |
| Alarm | High temperature alarm, visual (red LED) and acoustic (continuous buzzer) |

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| Nominal voltage [V] | Current consumption [A] | | Fuses IEC 127-2/V |
|---------------------|----------------------------|----------|----------------------|
| | at 50 Hz | at 60 Hz | |
| 230 | 2.3 | 2.4 | T4H250V |
| Other versions | | | |
| 100 | 5.9 | 5.9 | T8H250V |
| 110 | 4.9 | 5.2 | T8H250V |
| 115 | 4.7 | 5.0 | T8H250V |
| 127 | 4.2 | 4.5 | T8H250V |

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Electromagnetic
compatibility EMC

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II a

Classification
as per EC Directive 93/42/EEC
Annex IX

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10-972

UMDNS-Code
Universal Medical Device
Nomenclature System -
Nomenclature for medical
products

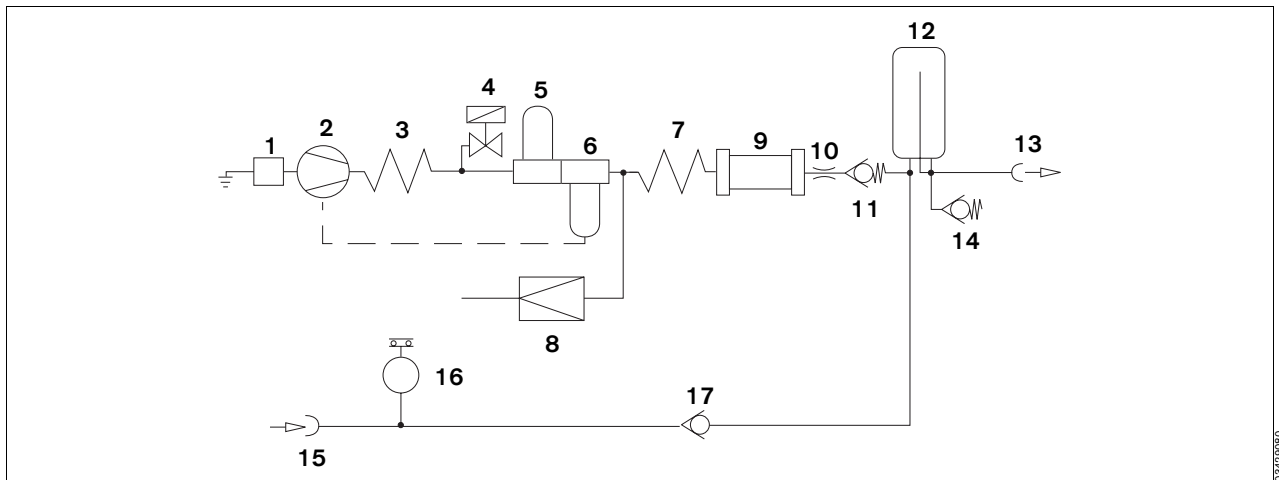
10-972

¹⁾ NIST = Non Interchangeable Screw Thread
(unvertauschbarer Anschluss)
²⁾ DISS = Diameter Indexed Safety System
(unvertauschbarer Anschluss)

¹⁾ NIST = Non Interchangeable Screw Thread
²⁾ DISS = Diameter Indexed Safety System

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- 1 Intake filter
- 2 Compressor unit
- 3 Cooling coil
- 4 Solenoid valve
- 5 Prefilter
- 6 Condensate trap
- 7 Heat exchanger
- 8 Pressure limiter
- 9 Diaphragm drier
- 10 Restrictor nozzle
- 11 Non-return valve
- 12 Pressure vessel
- 13 Medical air outlet
- 14 Relief valve
- 15 Input for standby mode
- 16 Pressure switch
- 17 Non-return valve

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Ambient air is drawn in via the intake filter 1, compressed in the compressor unit 2 and cooled in the cooling coil 3. The system is briefly opened by solenoid valve 4 when the compressor starts up.

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The compressed air is cleaned by prefilter 5. Condensate is collected and removed in the condensate trap 6.

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The purified air is re-heated by the heat exchanger 7 in order to prevent further condensation forming before it reaches the diaphragm drier 9.

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In the diaphragm drier 9, the air is dehumidified to a dew point at least 5 °C below the ambient temperature.

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The dehumidified air flows via the restrictor nozzle **10** and non-return valve **11** to the pressure vessel **12**, where it is withdrawn via the coupling **14**.

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In standby mode, the ventilator draws air from the central medical air supply via the coupling **15**, non-return valve **17** and self-closing coupling **13**. The compressor unit **2** is on standby.

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If the pressure in the central medical air supply drops below $2.7-0.3 \text{ bar}^{1)}$, the compressor unit **2** is activated via pressure switch **16**.

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When the pressure in the central medical air supply reaches $3.4 \text{ bar}^{2)}$ again, the compressor unit **2** is switched off again by the pressure switch **16**.

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Pressure regulator **8** limits the system pressure to 4 bar and relief valve **14** protects the unit against excessively high pressure from the central medical air supply.

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Order List

| Description | Order No. |
|---|-----------|
| Medical air compressor system | 84 14 350 |
| Medical air compressor | |
| Standard version 230 V, 50 Hz | 84 13 890 |
| Option for standby mode | 84 13 939 |
| Pressure hoses for operating the compressor on the Evita trolley: | |
| NIST connecting hose 0.8 m, black | 84 11 517 |
| NIST connecting hose 0.8 m, gas colour-coding to ISO 32 | 84 11 538 |
| DISS connecting hose 0.8 m, yellow | 84 14 454 |
| Set of replacement parts for the 6000 hour service for medical air compressors from ARMJ 0020 (see rating plate) | 84 14 501 |
| The following hoses may be ordered instead of the 0.8 m hose: | |
| NIST connecting hose 5 m, black | 84 11 519 |
| NIST connecting hose 5 m, gas colour-coding to ISO 32 | 84 11 541 |
| For standby mode: | |
| NIST medical air connecting hose 3 m | M 34 408 |
| NIST medical air connecting hose 5 m | M 34 409 |
| Trolley | 84 11 520 |
| Technical documentation on request | |

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© Dräger Medical AG & Co. KG